

## Foreword

Thanks you for choosing the foot massager to complement your healthcare needs. This product has been tested to be safe and reliable for personal use. Nevertheless, please read this user manual carefully in order to enjoy the product benefits fully.


## Safety Precautions

### Forbidden: ⚠ Warning

( It may cause injures and deaths if against the belowing rules. )



Please do not allow children play on the machine, or may cause injured or trouble.

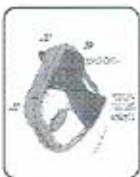

Please do not use it in the wet or high temperature area (especially saunas and bathroom), or may cause electric shock, fire disaster may happened or product become mildewed.



Please consult the doctor if the children want to use , or the electric shock or injured may happened.

1. Please place the device horizontally to ensure it works safely.
2. Please take out the adapter if unexpected power outages.
3. The auto timing is 15 minutes, 15-30 minutes massage each day can release the muscle ache or nerve ache.
4. Please do not scratch the massager cover by a knife or sharp tools.




## Specification

Name: Shiatsu Foot Massager	Model: AS-617	Safety design: <input type="checkbox"/>
Voltage: DC12V	Power: 24W	State of Origin: China
G.W.: 4.0kg	N.W.: 3.5kg	Production Date:
Measurement: 371x369x179mm	QTY: 1pcs	

## Operation and functions

1. Press the Power button. The indicator shining, the machine is ready for working.
2. Press the Auto button, it starts 3 modes automatically.
3. Press the Intensity, Speed, Direction button manually to choose different function separately.
4. Press the Power Button again to stop the massage.
5. Auto Timing is 15 minutes.

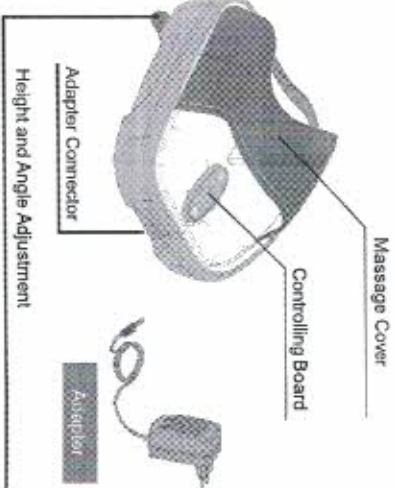

**Message Intensity**  
Low/Middle/High  
Intensity is optional

**Speed**  
Slow/Middle/Fast  
speed is optional

**ON/OFF**

**Auto Mode**  
3 auto mode  
adjusting function

**Heating**  
Turn on or off the  
heat function

**Direction**  
Clockwise and  
anticlockwise optional



## Care & Maintenance

1. Spot clean with a damp cloth.
2. DO NOT machine wash or dry clean.
3. DO NOT clean while switched on.
4. DO NOT leave in full sunlight or in a damp place.
5. Always switch power to off position after use.
6. After cleaning, store unit in a dry and cool place.